

RPA547Hu01 50µg
Vascular Cell Adhesion Molecule 1 (VCAM1)
Organism Species: Homo sapiens (Human)
Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

8th Edition (Revised in Jun, 2013)

[PROPERTIES]

Residues: Gln109~Ser318 (Accession # P19320),
with N-terminal His-Tag.

Host: *E. coli*

Subcellular Location: Membrane; Single-pass
type I membrane protein.

Purity: >95%

Endotoxin Level: <1.0EU per 1µg
(determined by the LAL method).

Formulation: Supplied as lyophilized form in PBS,
pH7.4, containing 5% sucrose, 0.01% sarcosyl.

Predicted isoelectric point: 5.5

Predicted Molecular Mass: 24.8kDa


Accurate Molecular Mass: 27kDa as determined by SDS-PAGE reducing conditions.

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

Note: The possible reasons that the actual band size differs from the predicted are as follows:

1. Splice variants: Alternative splicing may create different sized proteins from the same gene.
2. Relative charge: The composition of amino acids may affects the charge of the protein.
3. Post-translational modification: Phosphorylation, glycosylation, methylation etc.
4. Post-translation cleavage: Many proteins are synthesized as pro-proteins, and then cleaved to give the active form.
5. Polymerization of the target protein: Dimerization, multimerization etc.

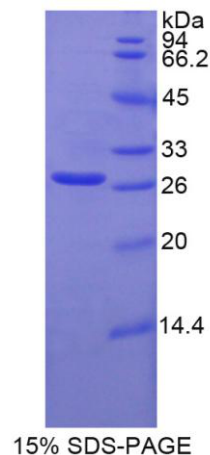

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The target protein is fused with N-terminal His-Tag, its sequence is listed below.
MGHHHHHHS~~G~~SEF-QV EIYSFPKDPE IHLSGPLEAG KPITVKCSVA DVYPFDRLEI
DLLKGDHLMK SQEFLEDADR KSLETKSLEV TFTPVIEDIG KVLVCRALKH IDEMDSVPTV
RQAVKELQVY ISPKNTVISV NPSTKLQEGG SVTMTCSSEG LPAPEIFWSK KLDNGNLQHL
SGNATLTLIA MRMEDSGIYV CEGVNLIGKN RKEVELIVQE KPFTVEIS

[REFERENCES]

1. Kim K.H., *et al.* (2012) *Phytother Res* 26:1037-1047.
2. Kanmaz-Ozer M., *et al.* (2012) *Eur. J. Obstet. Gynecol. Reprod. Biol.* 160:174-178.
3. Cartee T.V., *et al.* (2012) *J. Dermatol. Sci.* 65:86-94.
4. Saikia U.N., *et al.* (2012) *Cardiovasc. Pathol.* 21:83-88.



Gentaur Molecular Products
Voortstraat 49
1910 Kampenhout, BELGIUM
Tel 0032 16 58 90 45 | Fax 0032 16 50 90 45
www.gentaur-worldwide.com
info@gentaur.com